



+ Special complications

Model	Effective Length of the Needle Tube	Length of the Localization String
SS510-10	100mm	100mm
SS510-10L	100mm	135mm
SS510-12	120mm	100mm
SS510-12L	120mm	155mm
SS510-15	150mm	100mm
SS510-15L	150mm	185mm

+ Product Features

Precise Localization

Precise puncture and location under CT guidance. Multiple claws form gives the localization hook the best reliable anchoring characteristics in the lesion sites.

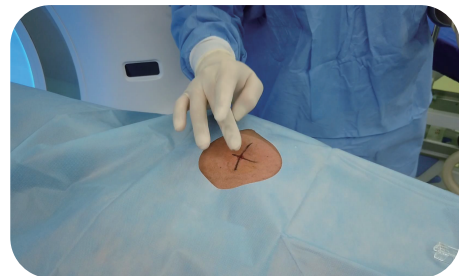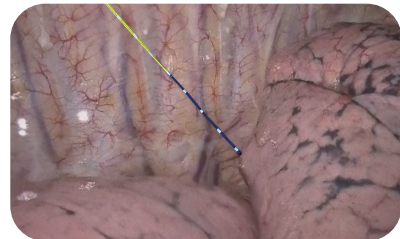

Anchor Hook Design for Secure Fixation

Features a four-prong rounded barbed configuration for reliable tissue engagement, with blunt tips that effectively reduce the risk of bleeding.

Flexible Thread Localization for Low Migration Risk

The flexible thread effectively buffers traction forces, minimizing anchor migration and ensuring more stable positioning.

Visual Depth Indication for Intuitive Identification

Color-coded depth markings allow for easy identification of nodule depth, facilitating accurate intraoperative judgment.

Disposable Pulmonary Nodules Localization Needle

Senscure 胜杰康
人杰质胜 健企康民

